

MINI POCKET BEAUTICIAN

THE SKIN SCRUBBER



USER'S MANUAL

HS-CO20

SMALLEST SIZE AMONG THE MARKET. DESIGNED ESPECIALLY FOR NOSE PART.

Function

Vibration: 23,000 ~ 25,000 times per second to turn the water drop into tiny particles to cleanse the blackheads on nose part.

Working Theory


We use PZT ceramics to turn the electricity into vibration power. It atomizes the water drop to cleanse dirt inside the nose skin pores.

Use Instruction

Spray water or cleansing liquid averagely on nose/facial skin.

Use the Mini Pocket Beautician with 30~45 degrees to slightly contact the skin to achieve deep cleansing.

It not only cleanses blackheads, but also helps remove the makeup residual, acnes and aged skin.

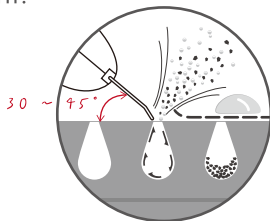


Attention

Read The Following Instructions Carefully

1. Please do not use the device in the bathroom or other humid environment.
2. Please do not crash or twist the device while using.
3. Please do not stay in the same position too long while using.
4. Do not press the spatula hard into the skin.
5. Use the earplug if the sound from the spatula makes the user feel uncomfortable.
6. Keep the spatula in the dry and cool place and do not wash the device directly under water.
7. Ask the doctor if you have any skin problem before using.
8. Do not use the device while the skin is complete dry. It is suggested to use when the skin is moistured..

Button & Indicator Light Description

1. Power : On / Off the indicator is in white color.
 - One Press to turn on the device. The indicator lights up in white color.
Press again the indicator light disappears and the device turns off.
One Button, One touch to operate.
2. Operation Time : 3 minutes each round.

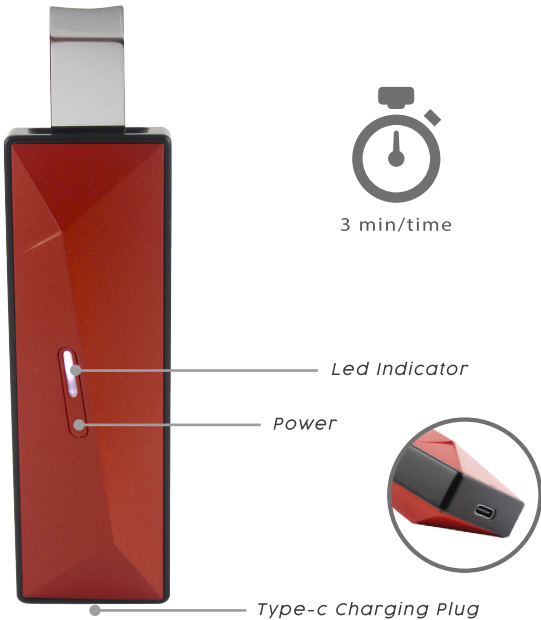

Steps to Use

Use slightly and no need to have much strength

■ Deep Cleansing

1. Cleanse the skin.
2. Spray the water on your nose/ face.
3. Use the skin scrubber with the angle of 30 to 45 degrees compared to our skin.
4. Keep not much water on skin scrubber to insure the cleansing effect.
5. Do not use the skin scrubber if your skin has no water on it.

Lead-in Mode

1. Calm down the skin after cleansing
2. Use the other side of the skin scrubber to message the nose/face skin.
3. Use along with the skin care products.
4. Pay attention to the cleansing direction and angle.

Power & Charging Indication

- Use the Type-C plug with cell phone power supply. The output voltage of the power has to be 5V.
- Charing Indication
 - Charging :
The power indicator would flash slowly to show the charging status.
 - Charging Completion :
The indicator light stops flashing and no light.
- Use Time : 30 minutes/ one fully charging
- Charging Time : 2.5 hrs.

Specification

Model : C020

Input/Output : AC100-240V / DC 5V , 0.8A

Size : L 137 x W 35 x H 30 mm

Weight : 55 g

Frequency : 23,000~25,000 Hz